
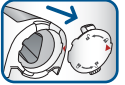



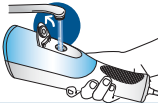


QUICK START GUIDE

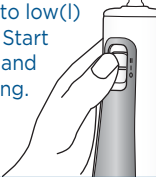
Review full instruction book prior to first use

1 FIRST USE


1. Remove reservoir. 
2. To remove battery door, twist, and pull. 
3. Insert 3AA batteries. Included 
4. To replace battery door, line-up arrows, push and twist. 
5. Insert tip firmly until it clicks into place. 
6. Fill reservoir with warm water. 

2 START WATER FLOSSING

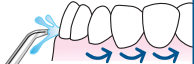
TO AVOID MESS: PLACE THE TIP IN YOUR MOUTH BEFORE TURNING UNIT ON

Slide power switch to low (I) or high (II) position. Start with lowest setting and increase to your liking. 

3 LEAN OVER SINK

Close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink. 

4 TECHNIQUE

Aim water at the gumline at a 90° angle. Follow the gumline and pause briefly between teeth. 

5 HOW TO VIDEO




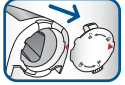


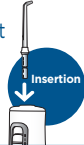

Scan here or visit waterpik.com/OH03

Register your product and receive over \$30 in savings. waterpik.com/registration

GUIDE DE DÉMARRAGE RAPIDE

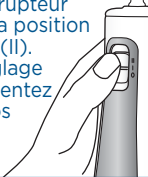
Lisez l'intégralité du manuel d'instructions avant utilisation

1 PREMIÈRE UTILISATION


1. Retirez le réservoir. 
2. Pour retirer le couvercle du compartiment des piles, tournez et tirez. 
3. Insérez 3 piles de type AA. Incluses 
4. Pour réinstaller le couvercle du compartiment des piles, alignez les flèches, poussez et tournez. 
5. Insérez fermement l'embout jusqu'à ce qu'il s'emboîte. 
6. Remplissez le réservoir d'eau tiède. 

2 DÉMARRAGE DE L'HYDROPROPULSEUR


POUR ÉVITER LES ÉCLABOUSSURES : PLACEZ TOUJOURS L'EMBOUT DANS LA BOUCHE AVANT DE METTRE L'APPAREIL EN MARCHÉ.

Faites glisser l'interrupteur d'alimentation sur la position basse (I) ou élevée (II). Démarrez sur le réglage le plus bas et augmentez la pression selon vos préférences. 

3 PENCHEZ-VOUS PAR-DESSUS L'ÉVIER

Gardez vos lèvres suffisamment fermées pour empêcher les éclaboussures, tout en laissant l'eau s'écouler de votre bouche dans le lavabo. 

4 TECHNIQUE

Dirigez le jet d'eau vers la ligne gingivale à un angle de 90°. Suivez la ligne gingivale et faites de courtes pauses entre les dents. 

5 VIDÉO D'INSTRUCTIONS



Scannez ce code ou visitez le site Internet waterpik.com/OH03

Enregistrez votre produit et recevez plus de 30 \$ d'économies. waterpik.com/registration