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Replace the reservoir onto the base. There should be no gap between the reservoir and the base.

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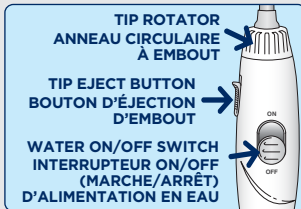
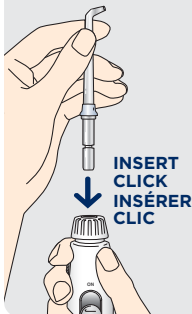
**Caution:** Do not slide tip eject button while unit is running.

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#### Prior to Very First Use of Unit:

Set pressure control dial to **highest position**, point the tip into sink. Turn power on. Slide water switch **ON** until water flows from the tip with substantial force. Turn power **OFF**.

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
For first time users, adjust pressure control to **lowest position**. Gradually increase pressure to your liking.

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Press the mode button once and one LED will display current mode. Press mode button again to change modes. Floss Mode is suggested for first time use.

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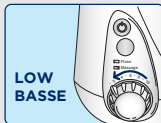
*continued on back / voir suite au dos*

This device is only compatible with a voltage system of 120VAC, 60Hz. Do not plug this device into a voltage system that is different.

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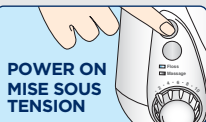
**POWER ON  
MISE SOUS  
TENSION**

**WATER ON  
EAU-MARCHE**

With water switch on handle in OFF position, turn power on. Lean over sink and close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink. Turn water switch ON. When finished, turn water switch and power OFF.

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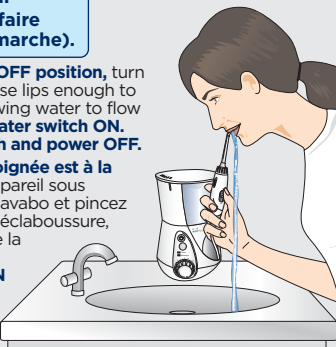During use, you can stop water flow by sliding water switch on handle to the OFF position.

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
Aim the water stream just above the gumline at a 90 degree angle.

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Starting with the back teeth, follow the gumline and pause briefly between teeth. Continue until all areas around and between teeth have been cleaned.

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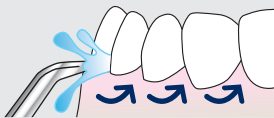**Timer:** In Floss Mode the timer will provide a short pause at 30 seconds for pacing and at 1-minute to indicate recommended cleaning time is complete.

**Note:** Timer may be deactivated if desired. See instruction book for details.

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**Massage Mode:** For best results, use the Floss Mode as instructed. Switch to Massage Mode, adjust pressure as desired, and aim the water stream at your gums and massage throughout your mouth.

**Note:** See instruction book for full details.

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